

Dornier
AXIS LEO

Clean & Clear, Every Time

Single-use digital flexible ureterorenoscope



9.6

Take
your
pick

7.5

AVAILABLE IN 2 SIZES



Dornier *MedTech*



If size matters,

make your choice without compromise on performance.

Available in 2 sizes,
7.5 Fr & 9.6 Fr

7.5 Fr

9.6 Fr


UAS

Ureteral Access Sheath



| Video Processor

- Multiple connectivity options (DVI, S-Video, CVBS) allow compatibility with most monitor types for easy integration into any OR
- User-friendly menu navigation
- Integrated USB data storage option

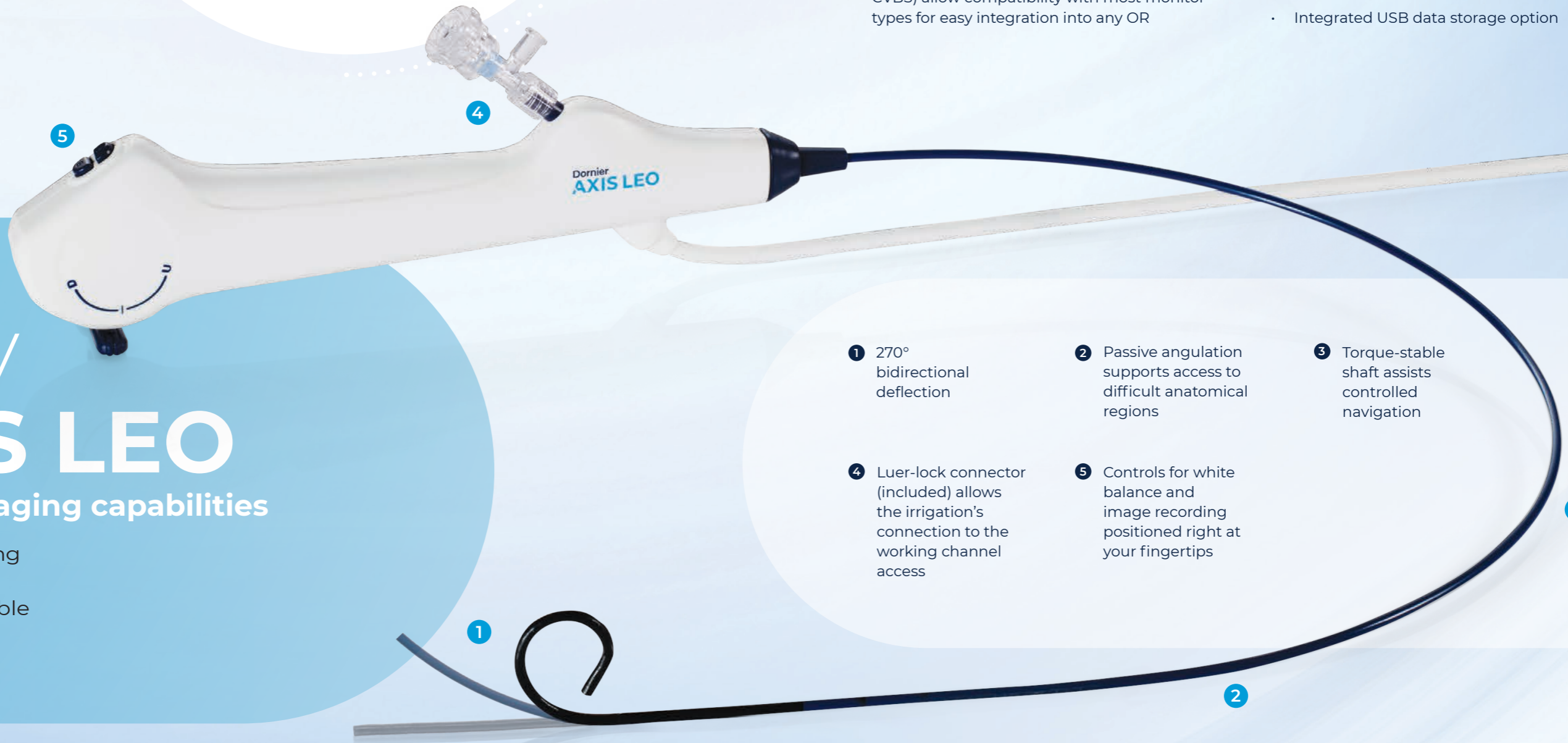

Why AXIS LEO

Brilliant imaging capabilities

- No reprocessing
- No repairs
- Readily available
- Always clean
- Cost-effective

- 1 270° bidirectional deflection
- 2 Passive angulation supports access to difficult anatomical regions
- 3 Torque-stable shaft assists controlled navigation
- 4 Luer-lock connector (included) allows the irrigation's connection to the working channel access
- 5 Controls for white balance and image recording positioned right at your fingertips

AXIS LEO Flexible Video Ureterorenoscope

	9.6 Fr	7.5 Fr
Order information	K1046369 9.6 Fr URS EU deflection	K1048308 7.5 Fr URS EU deflection
	K1046370 9.6 Fr URS US deflection	K1048309 7.5 Fr URS US deflection
Distal tip	9.6 Fr / 3.2 mm	7.5 Fr / 2.5 mm
Working channel	3.6 Fr / 1.2 mm	
Working length	650 mm	
Deflection	270° up / down	
Illumination	Optical Fiber	Dual LED
Field of view	120°	
Depth of field	3 - 50 mm	
CMOS pixel	400 * 400	

AXIS LEO Video Processor

Order information	K1046368
Dimensions	(W) 24 cm x (L) 30 cm x (H) 6 cm
Weight	0.75 kg
Connections	USB 2.0, DVI, S-Video, CVBS



Scan the QR code for our
global office locations

www.dornier.com/locations